



Hummingbird

Medical Workstation on Wheels



Laptop



Powered



Manual Lift



Amico

Patient



Improved Communication



Increased Treatment Compliance



Enhanced Security



Stronger Outcomes

Caregiver



Increased Workflow Efficiency



Fewer Mistakes



Less Pain



Reduced Fatigue

Technician



Maximal Uptime



Minimal Maintenance



Fleet Management



Multi-Level Support

Management



Reduced Liability




Better Patient Outcomes



Increased Caregiver Productivity



Superior Management and Integration

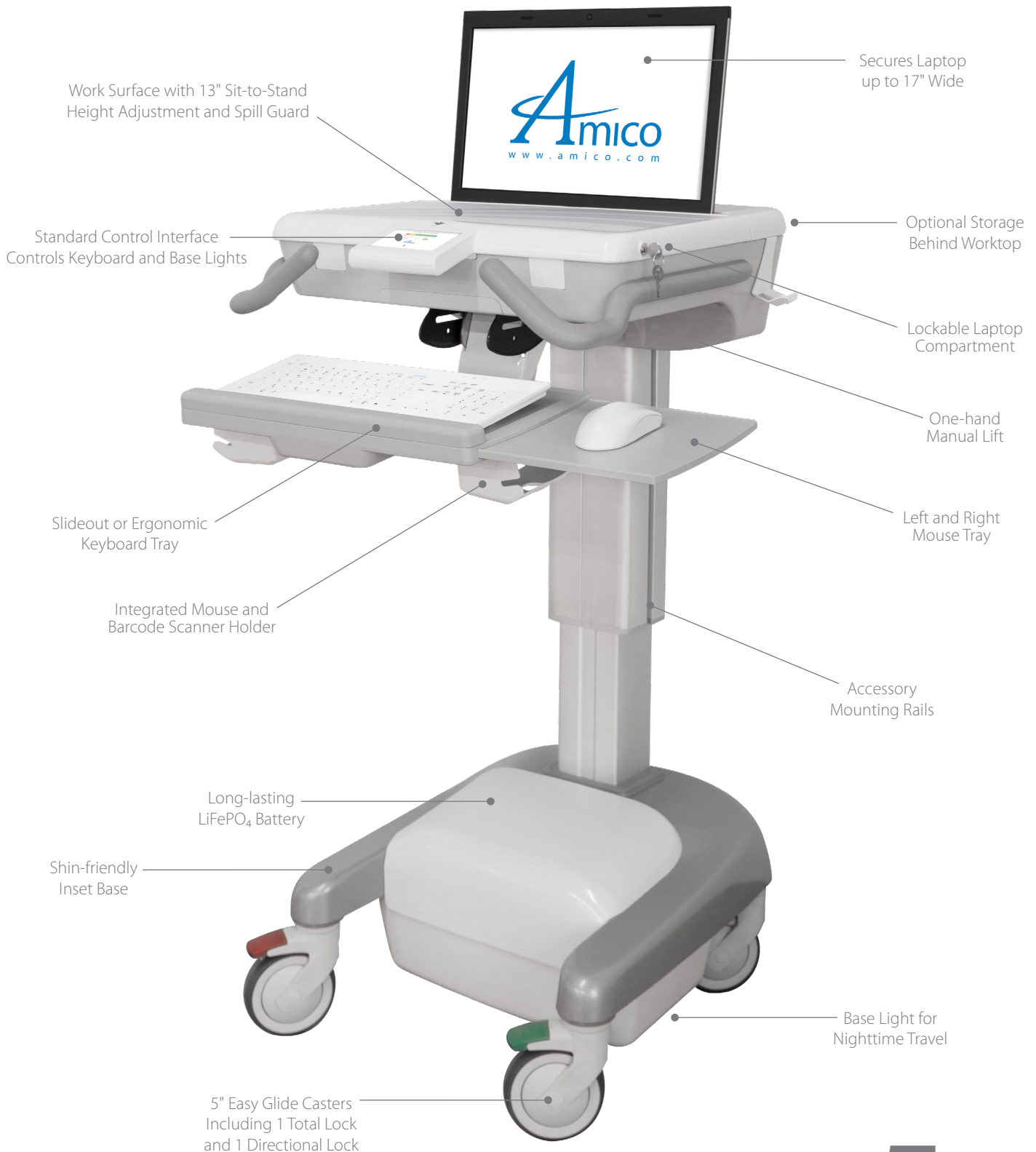

5YEAR
warranty

Hummingbird Standard Equipment



Worktop

18 sq" (1200 cm²) work surface with 8 fl. oz (236.59 ml) spill guard ideal for documentation and medication preparation



Standard Interface

Controls worktop height, cart lights, and provides battery status and alerts



Keyboard Light

Illuminated keyboard allows you to work at night without disturbing your patients



Ergonomic Handles

Optimally oriented handles make maneuvering your cart easy and reduces fatigue



Secure I.T.

Secure your CPU or Laptop in our lockable compartment




Barcode Scanner/Mouse Holder

Place your mouse or barcode scanner underneath the keyboard tray to prevent them from falling during travel




Power Cord

6' (1.83 m) coiled power cord for convenient connection to wall outlet

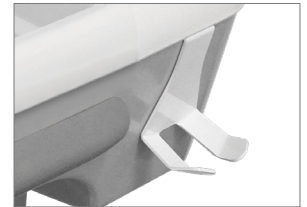


Cord Holder

Neatly store your power supply cord on the side of the cart

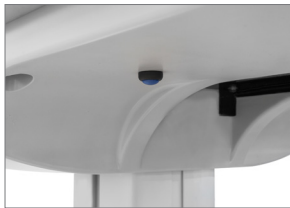

Power Reset

Integrated reset button for rebooting your computer system (requires BIOS setting)

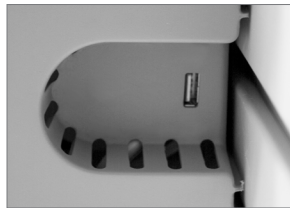


External USB Ports

2 hidden USB ports for connecting auxiliary devices



Accessory Rails

Integrated rails on both sides of the column to easily mount your storage accessories

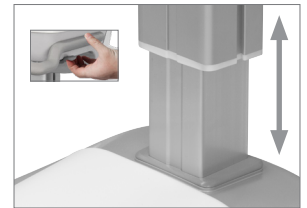


Manual Lift

Go from sitting to standing in seconds with our heavy-duty pneumatic lift: 13" (33.02 cm) worktop height adjustment range



Effortless Rolling

5" casters make travel effortless, including directional caster for enhanced steering and total lock caster to stop movement

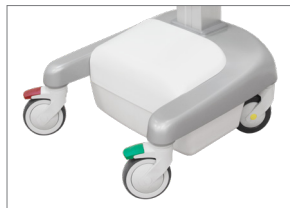

Legroom

U-Shape base saves your shins while traveling and gives users enough legroom while sitting



Base Light

See where you're going with an integrated LED light at the base of the cart



Plug and Play

Ships assembled: add your peripherals/accessories and your cart is ready to use



Fleet Management

Track the cart status of your entire fleet and receive critical service alerts in real time with BirdsEye Fleet Management Software

Hummingbird Customization Options



LiFePO₄ Battery

3-hour fast charge to ~90% full, 4 hours to full charge from empty
19.2Ahr (230.4WHr)
38.4Ahr (460.8WHr)
57.6Ahr (691.2WHr)




12' Power Cord

12' (3.66 m) extended length power cord for easy access to hard-to-reach outlets




Power System

Power connections available for all countries, choose AC-only or AC/DC power output

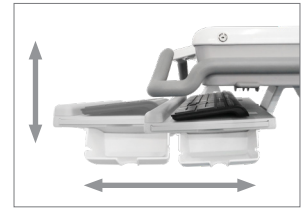

Slideout Keyboard

Pull out or tuck away your keyboard; includes right and left-handed mouse trays

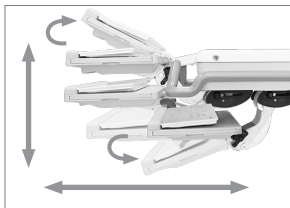

Ergonomic Keyboard Tray

4 ranges of motion such as; tilt, side-to-side, height, and in-and-out for optimal ergonomic positioning and patient orientation



Hummingbird Add-On Accessories

Bump Storage Accessories



Bump Storage Bar

The Bump Storage Rail mounts behind your worktop and lets you place frequently-used items close at hand



Universal Sharps Container*

Adjustable mount to accommodate common sharps containers



Hand Sanitizer Dispenser Mount*

Mounting plate for your sanitizer dispenser (ask sales rep about your dispenser type)



Equipment Wipes Mount*

Mounting bracket for your equipment wipes dispenser



Glove Box Holder*

Storage for boxes of gloves (works with tissue boxes too)


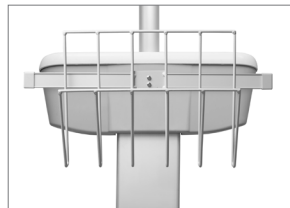


Chart Holder*

Locate your charts close at hand, behind the worktop

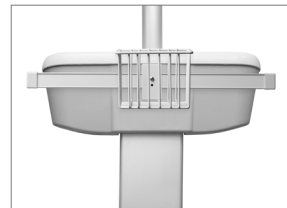

Wire Basket*

5" x 5" x 4" (12.7 cm x 12.7 cm x 10.16 cm) wire basket for storage



Diagnostic Station*

Add any combination of otoscope/ophthalmoscope station, specula dispenser, aneroid, and thermometer

Worktop Accessories



Worktop Light

Light up your worktop for nighttime charting without patient disruption



RFID

RFID embedded in the worktop facilitates secure PC access and enhanced infection control




USB Splitter

4-Port USB Splitter for extra devices



USB Extension

Extends the reach of your USB by 3' (91.44 cm)

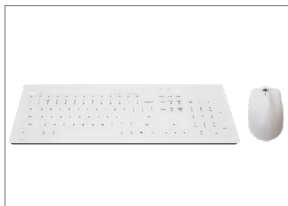

Antimicrobial Keyboard and Mouse

Antimicrobial keyboard and mouse (ask sales rep for options)



Custom Label

Personalize your carts with a worktop surface label of your choice



Built-in Ethernet Port

In work surface for network connections



Voltage Converter (110V to 220V)

Enable the powered Hummingbird Cart to operate on 220V / 230V power supply along with associated devices on the cart

Hummingbird Add-On Accessories

Vertical Rail Column Accessories



Storage Drawers

Key-locking drawer offers simple security for equipment and essential supplies



Baskets

Sterilizable Wire Baskets
10" x 6" x 9" (25.4 cm x 15.24 cm x 22.86 cm)
13" x 9" x 7" (33.02 cm x 22.86 cm x 17.78 cm)



Cuff Basket

Sterilizable Wire Baskets
13" x 9" x 7" (33.02 cm x 22.86 cm x 17.78 cm)



Wipes Holder

Mounting bracket for your equipment wipes dispenser




Printer Shelf

Open storage for your printers
small 6" x 6" (15.24 cm x 15.24 cm)
medium 12" x 12" (30.48 cm x 30.48 cm)
large 12" x 20" (30.48 cm x 50.8 cm)



Glove Box Holder

Storage for boxes of gloves
(works with tissue boxes too)

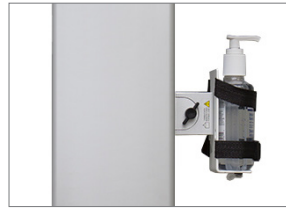

Sanitizer Holder

Securely store your sanitizer pump



Universal Sharps Holder

Adjustable mount to accommodate common sharps containers



Waste Bin

Store waste on the cart



Amico Accessories Inc. | 85 Fulton Way, Richmond Hill, ON L4B 2N4, Canada
Toll Free Tel: 1.877.264.2697 | Tel: 905.763.7778 | Fax: 905.763.8587
Email: info@amico-accessories.com | www.amico.com

Represented By: