



Hummingbird

Medical Workstation on Wheels



All-In-One,
LCD, or Tablet



Non-Powered



Manual Lift



Amico

Patient



Improved Communication



Increased Treatment Compliance



Enhanced Security



Stronger Outcomes

Practitioner



Increased Workflow Efficiency



Fewer Mistakes



Less Pain



Reduced Fatigue

Technician



Maximal Uptime



Minimal Maintenance



Fleet Management



Multi-Level Support

Management



Reduced Liability



Better Patient Outcomes




Increased Practitioner Productivity



Superior Management and Integration




Hummingbird

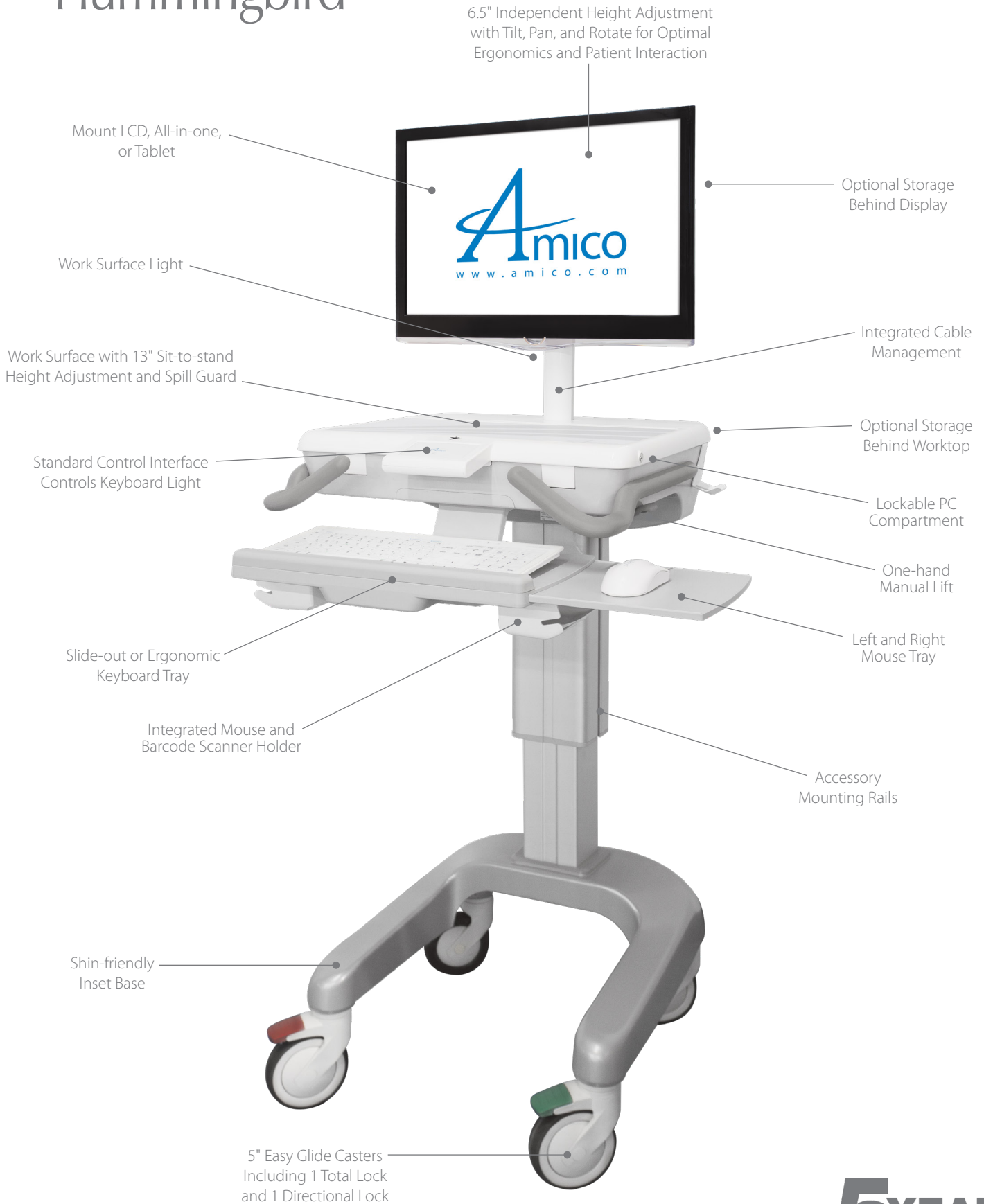


5YEAR
warranty

Hummingbird Standard Equipment

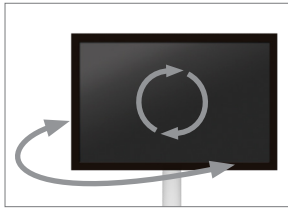

Display Movement

+/-90° screen panning and rotation for workflow, patient interaction, and visibility during travel

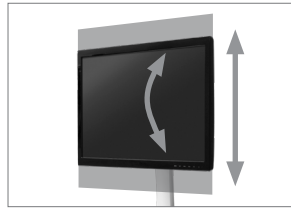

Display Ergonomic Adjustment

Independently adjust monitor height without tools for optimal ergonomic positioning: 6.5" (16.51 cm) height and +/- 15° tilt range



Cable Management

Conceal monitor and accessory cables inside the LCD post to keep them organized and clean



Worktop

18 sq" (1,200 cm²) work surface with 8 fl. oz (236.59 ml) spill guard ideal for documentation and medication preparation



Keyboard Light

Illuminated keyboard allows you to work at night without disturbing your patients



Ergonomic Handles

Optimally oriented handles make maneuvering your cart easy and reduces fatigue



Secure I.T.

Secure your CPU or Laptop in our lockable compartment



Barcode Scanner/Mouse Holder

Place your mouse or barcode scanner underneath the keyboard tray to prevent them from falling during travel



Power Cord

6' (1.83 m) coiled power cord for convenient connection to wall outlet



Cord Holder

Neatly store your power supply cord on the side of the cart




External USB Ports

2 hidden USB ports for connecting auxiliary devices



Accessory Rails

Integrated rails on both sides of the column to easily mount your storage accessories

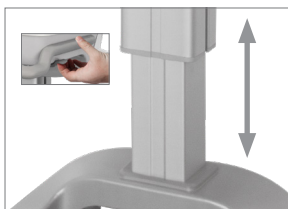


Manual Lift

Go from sitting to standing in seconds with our heavy-duty pneumatic lift: 13" (33.02 cm) worktop height adjustment range



Effortless Rolling

5" casters make travel effortless, including directional castor for enhanced steering and total lock castor to stop movement

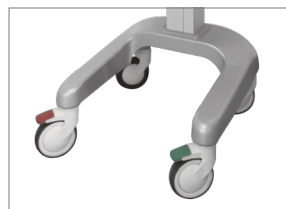

Legroom

U-shaped base saves your shins while traveling and gives users enough legroom while sitting



Plug and Play

Ships assembled: add your peripherals/accessories and your cart is ready to use


Hummingbird Customization Options




Standard Interface
Controls cart lights




Power System
Power connections available for all countries

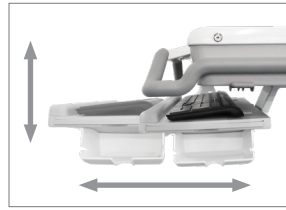

Slide-out Keyboard
Pull out or tuck away your keyboard; includes right and left-handed mouse trays

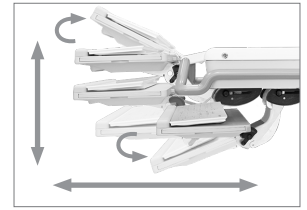

Ergonomic Keyboard Tray
4 ranges of motion such as; tilt, side-to-side, height, and in-and-out for optimal ergonomic positioning and patient orientation

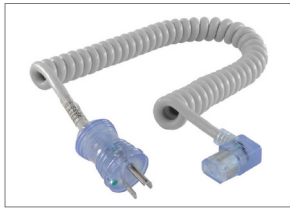

12' Power Cord
12' (3.66 m) extended length power cord for easy access to hard-to-reach outlets



Display (LCD) Mount
VESA Mount for your display




AIO Mount
VESA Mount for your all-in-one



Tablet Mount
VESA Mount for your tablet enclosure

Hummingbird Add-On Accessories

Vesa Accessories

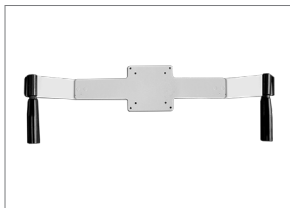

Display Handles
Ergonomic handles make screen adjustment more comfortable for practitioners



CPU Mount
Mount your small sized CPU behind the display



Tablet Enclosure
Please contact your sales representative to discuss a range of enclosure options

Hummingbird Add-On Accessories

Perch Storage Bar Accessories




Perch Bar

The Perch Storage Rail mounts behind your display and lets you place frequently-used items close at hand



Barcode Scanner

Mount your wired or wireless barcode scanner beside your display

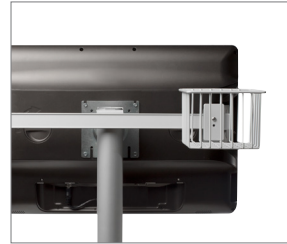


5" x 5" x 4" Basket

5" x 5" x 4" (12.7 cm x 12.7 cm x 10.16 cm) wire basket for storage



Sharps Container*

Adjustable mount to accommodate common sharps containers

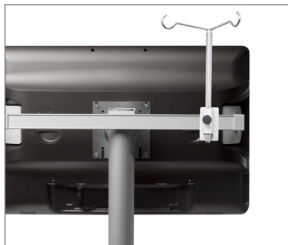

I.V. Hook

Hang I.V. bags behind your screen (2 or 4 hook options available)



Hand Sanitizer Pump Holder*

Securely store your hand sanitizer pump



Hand Sanitizer Dispenser Mount*

Plate for mounting your sanitizer dispenser (ask sales rep about your dispenser type)



Equipment Wipes Holder*

Mounting bracket for your equipment wipes dispenser



Printer Mount*

Mounting plate accommodates commonly used small printers including Zebra QL220

*Contact your sales rep for options and availability

Hummingbird Add-on Accessories

Bump Storage Accessories



Bump Storage Bar

The Bump Storage Rail mounts behind your worktop and lets you place frequently-used items close at hand



Universal Sharps Container*

Adjustable mount to accommodate common sharps containers



Hand Sanitizer Dispenser Mount*

Mounting plate for your sanitizer dispenser (ask sales rep about your dispenser type)



Equipment Wipes Mount*

Mounting bracket for your equipment wipes dispenser



Glove Box Holder*

Storage for boxes of gloves (works with tissue boxes too)


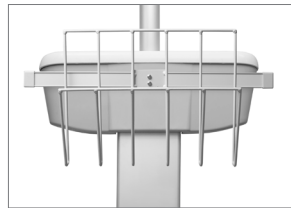


Chart Holder*

Locate your charts close at hand, behind the worktop

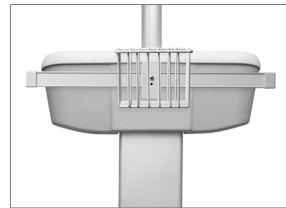

Wire Basket*

5" x 5" x 4" (12.7 cm x 12.7 cm x 10.16 cm) wire basket for storage



Diagnostic Station*

Add any combination of otoscope/ophthalmoscope station, specula dispenser, aneroid, and thermometer

Worktop Accessories



Worktop Light

Light up your worktop for nighttime charting without patient disruption



RFID

RFID embedded in the worktop facilitates secure PC access and enhanced infection control




USB Splitter

4-port USB Splitter for extra devices



USB Extension

Extends the reach of your USB by 3' (91.44 cm)

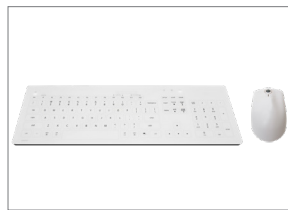

Antimicrobial Keyboard and Mouse

Antimicrobial keyboard and mouse (ask sales rep for options)



Custom Label

Personalize your carts with a worktop surface label of your choice



Built-in Ethernet Port

In work surface for network connections

*Contact your sales rep for options and availability

Hummingbird Add-On Accessories

Vertical Rail Column Accessories



Storage Drawers

Key-locking drawer offers simple security for equipment and essential supplies



Baskets

Sterilizable Wire Baskets
10" x 6" x 9" (25.4 cm x 15.24 cm x 22.86 cm)
13" x 9" x 7" (33.02 cm x 22.86 cm x 17.78 cm)



Cuff Basket

Sterilizable Wire Baskets
13" x 9" x 7" (33.02 cm x 22.86 cm x 17.78 cm)



Wipes Holder

Mounting bracket for your equipment wipes dispenser



Printer Shelf

Open storage for your printers
Small 6" x 6" (15.24 cm x 15.24 cm)
Medium 12" x 12" (30.48 cm x 30.48 cm)
Large 12" x 20" (30.48 cm x 50.8 cm)



Glove Box Holder

Storage for boxes of gloves
(works with tissue boxes too)



Sanitizer Holder

Securely store your sanitizer pump



Universal Sharps Holder

Adjustable mount to accommodate common sharps containers



Waste Bin

Store waste on the cart