

Horizontal Bed Head Trunking Unit with Lights

General Specifications

Fully customized to meet your facility's specific needs.

Designed in accordance with NFPA and/or HTM design guides for Bed Head units.

Powder coated access panels and laminate colors for product longevity and to meet your room décor.

Equipped with LED lights with a 10 year warranty.

Seamless overlapping edges.

Optional accessory rails.

Product Applications


Enclosure shall be constructed of extruded, anodized aluminum alloy sections to provide a modular surface mounted unit.

The Sapphire Series LT system can be manufactured in lengths that range from 1,220 mm to 2,440 mm.

Fascia panels, top and bottom cover panels shall be removable to provide access for easy installation and maintenance of the services.

Laminated Fascia

Aluminum strips with High Pressure Laminate (HPL) with three (3) colors to choose from.

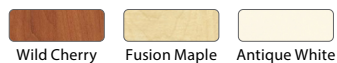

Medical Gas Outlets

Amico Corporation manufactures several compatible connection types for Medical Gas and Vacuum Outlets.

Electrical Wiring

Wire for standard and critical branch power circuits shall be #10 or #12 (as specified) type THHN stranded copper wire, 600 volt, with heat resistant thermoplastic insulation for hot (black) and neutral (white). Grounds shall be #10 type R THHN stranded copper wire (green). All ground conductors shall be factory installed in conduit.

Grounding and Bonding

Each power receptacle shall have a ground conductor connected to a grounding screw. A grounding bus will be installed in each device location.

Optional Accessory Rail



Features


- Any NFPA or HTM outlets can be used
- Removable laminate panels for access to electrical services
- Accommodates most configurations of medical gas outlets, electrical and communication services
- LED factory installed lights located at the top and bottom of the unit
- Three standard colors of hospital grade laminate to choose from
- Surface mounted medical gas outlets allow ease of access for annual maintenance
- Individual gas outlet connection
- Individual electrical wiring connection

Submittals and Verifications

Job specific shop drawings shall be produced for each project. As-Built drawings will be produced and provided at the time of shipment.

Amico Corporation will manufacture the equipment as per signed approval drawings and verification documents, which are provided in the job specific submittal package.

Horizontal Bed Head Trunking Unit with Lights

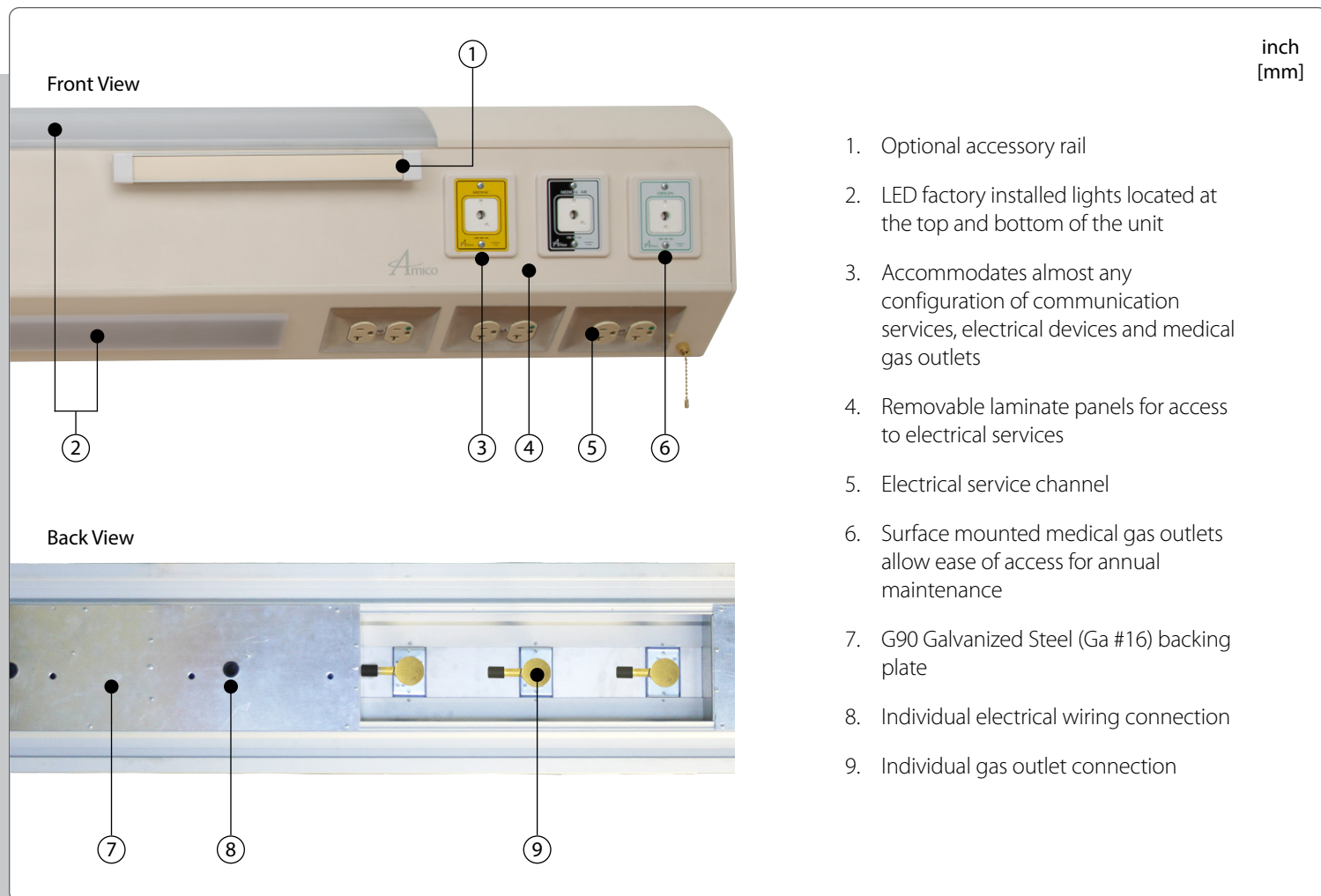

inch
[mm]

1. Optional accessory rail
2. LED factory installed lights located at the top and bottom of the unit
3. Accommodates almost any configuration of communication services, electrical devices and medical gas outlets
4. Removable laminate panels for access to electrical services
5. Electrical service channel
6. Surface mounted medical gas outlets allow ease of access for annual maintenance
7. G90 Galvanized Steel (Ga #16) backing plate
8. Individual electrical wiring connection
9. Individual gas outlet connection

The Letters "XX" Represent the Length of the Console:

- 48 = 1,220 mm*
- 60 = 1,525 mm
- 72 = 1,830 mm
- 84 = 2,140 mm
- 96 = 2,440 mm

*1,220 mm consoles can only house one (1) light, either on the top of the unit or the bottom.

Model Numbers

HWXX-09P-X-SGEN

09P = Sapphire Series Plus

The Letter "X" Represents the Number of Lights:

- L = 2 Lights (top and bottom)
- S = 1 Light (top or bottom)
- D = LED Lights (top and bottom)
- E = LED Light (top or bottom)
- 0 = No Lights

Distributed By: