

iCE Series Saturn Boom




Features

- Directly connected to the ceiling.
- Choice of fixed or 340° rotating option.
- Multiple console sizes available.
- Customizable outlet selection includes choice of electrical outlets, medgas services and data outlets.
- Hinged back door provides easy access to all internal components allowing for hassle-free maintenance.
- Integrated rails located on all corners of console allow for mounting of multiple accessories.
- Topcap prevents dust and bacteria from collecting on the top surface of the console with a round, easy to clean cover.

The iCE Series Saturn Booms offer a fixed alternative to the iCE Series. This unique series is directly connected to the ceiling and provides a fixed or rotating option (340°) to the user. The iCE Series Saturn Boom is completely customizable including everything from console size, colour, gases, outlets and accessories depending on your unique clinical needs. iCE Series Saturn Booms are available from 20 inch to 60 inch sizes.

Technical Specifications

Optional Rotation

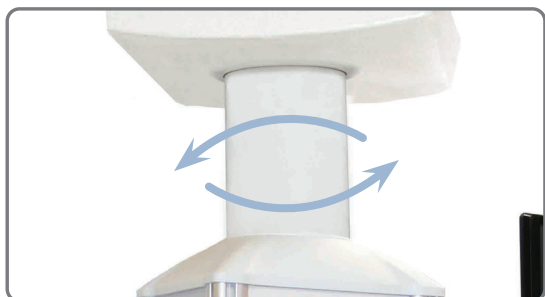

Integrated Equipment Rails



Direct Mount AHM Monitor Arm



General Specifications:

Input Voltage: 100 - 240 VAC

Power Consumption:

Ceiling and Adapter Bearing: 200 W

Console Bearing: 80 W

Switch Operation: Momentary

Braking System: Brake switches are provided on the ergonomic control handle on the console along with simplified operation diagram

Brake Type: Spring actuated friction brake

Torque: 520 lbs-in

Materials:

Arms: High Strength Powdercoated Aluminum

Bearing Housings: High Strength Steel

Console: Aluminum casing with Steel panels